# Prospektüs Prospectus

## POLISIL

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NONABSORBABLE, SURGICAL SUTURE, U.S.P./ Eur.Ph./TS-55303

**DESCRIPTION**

POLISIL polyester suture is a nonabsorbable, braided, sterile, surgical suture composed of Poly (ethylene terephthalate). It is prepared from fibers of high molecular weight, long-chain, linear polyester. POLISIL suture is uniformly coated with silicone elastomer or uncoated. The highly adherent coating is a relatively nonreactive. Coating compound which is nonabsorbable, acts as a lubricant and mechanically improves the physical properties of the braided uncoated suture by improving easy passage through tissues and overall handling qualities. POLISIL suture is braided for optimal handling properties and dyed green to increase visibility. The dye is D&C Green #6. POLISIL suture meets all requirements for nonabsorbable surgical sutures established by the United States Pharmacopoeia (U.S.P) and European Pharmacopoeia (E.P.).

**INDICATIONS**

POLISIL suture is indicated for use in general soft tissue approximation and/or ligation, including use in cardiovascular, ophthalmic and neurosurgical procedures.

**ACTIONS**

POLISIL suture elicits a minimal acute inflammatory reaction in tissue, followed by gradual encapsulation of the suture by fibrous connective tissue. Both polyester fiber suture material and the silicone coating are pharmacologically inactive.

**CONTRAINDICATIONS**

Not known.

**WARNINGS**

Users should be familiar with surgical procedures and techniques involving non-absorbable sutures before employing POLISIL suture for wound closure, as risk of wound dehiscence may vary with the site of application and the suture material used.

Do not re-sterilize. Discard open and unused sutures.

As with any foreign body, prolonged contact of any suture with salt solutions, such as those found in the urinary or biliary tracts, may result in calculus formation.

Acceptable surgical practices must be followed for the management of contaminated or infected wounds.

**PRECAUTIONS**

In handling POLISIL or any other suture material, care should be taken to avoid damage to needle and suture. Avoid crushing or crimping damage due to application of surgical instruments such as forceps or needle holders.

Multicurl sutures, which are treated to enhance handling characteristics requires the accepted surgical technique of flat and square ties with additional throws as warranted by surgical circumstances and the experience of the surgeon.

To avoid damaging needle points and swage areas, grasp the needle in an area one-third (1/3) to one-half (1/2) of the distance from the swaged end to the point. Reshaping needles may cause them to lose strength and be less resistant to bending and breaking. Users should exercise caution when handling surgical needles to avoid inadvertent needle sticks.

Dispose of contaminated and unused products are in accordance with local and facility requirements. Discard used needles in "sharps" containers.

**ADVERSE REACTIONS**

Adverse effects associated with the use of this device include wound dehiscence, calculi formation in urinary and biliary tracts when prolonged contact with salt solutions such as urine and bile occurs, infection, minimal acute inflammatory tissue reaction, and transitory local irritation at the wound site.

**HOW SUPPLIED**

POLISIL sutures are available as sterile, braided, green, and undyed (white) strands in U.S.P. sizes 7/0 to 5 (metric sizes 0.5 - 7) in a variety of lengths, with and without permanently attached needles. POLISIL sutures, green and undyed, are also available attached to Polytetrafluoroethylene (PTFE) pledget. POLISIL sutures are available in half, one, two, and three dozen boxes. POLISIL suture is for single use only.

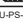
**STORAGE**

Store below 25 °C and keep away from sunlight.

Protect from humidity.

Do not use after expiry date.

**SYMBOLS USED ON LABELLING**

	Do not reuse!	<b>REF</b>	Catalogue Number
	Do not re-sterilize		Store below 25°C
	Do not use if pack age is damaged		Keep away from sunlight
	Manufacturer		Protect from humidity
	YYYY MM Date of Manufacture, Year		Recyclable pack
	YYYY-MM Expiry Date, Year - Month		Attention, See instruction for use
<b>STERILE EO</b>	Sterile EO: Ethylene oxide		<b>1014</b>
<b>LOT</b>	Batch Number		
	Dyed, Non-absorbable, Braided, Coated		

EU-PS-rev-03-05-05-2017

Issue date: 11.09.2012

"EASSI (The European Association of Surgical Suture Industry) has developed a system of symbols which is designed to describe various suture product characteristics in an intuitive, pictorial manner. The use of symbols is permitted by the Medical Device Directive (MDD 93/42/EEC) and enables companies to provide information to the customer without having to provide multilingual translations."

**BOZ TIBBI MALZEME SANAYI VE TICARET A.Ş.**  
Sağlık Sokak No:33/5 Sıhhiye/ANKARA/TÜRKİYE  
Tel:+90 (312) 435 50 80-81 Faks:+90 (312) 435 41 26  
web: www.boztibbi.com e-mail: boz@boztibbi.com

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